ABclonal www.abclonal.com

COX2/PTGS2 Rabbit pAb

Catalog No.: A1253 142 Publications

Basic Information

Observed MW

72kDa

Calculated MW

69kDa

Category

Polyclonal Antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human, Mouse, Rat

Background

Prostaglandin-endoperoxide synthase (PTGS), also known as cyclooxygenase, is the key enzyme in prostaglandin biosynthesis, and acts both as a dioxygenase and as a peroxidase. There are two isozymes of PTGS: a constitutive PTGS1 and an inducible PTGS2, which differ in their regulation of expression and tissue distribution. This gene encodes the inducible isozyme. It is regulated by specific stimulatory events, suggesting that it is responsible for the prostanoid biosynthesis involved in inflammation and mitogenesis.

Recommended Dilutions

WB 1:500 - 1:1000

IHC-P 1:50 - 1:200

Recommended starting **ELISA**

> concentration is 1 µg/mL. Please optimize the concentration based on your specific

assay requirements.

Immunogen Information

Gene ID Swiss Prot 5743 P35354

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 375-604 of COX2/PTGS2 (NP_000954.1).

Synonyms

COX2; COX-2; PHS-2; PGG/HS; PGHS-2; hCox-2; GRIPGHS; COX2/PTGS2

Contact

www.abclonal.com €


Product Information

Source **Isotype Purification** Rabbit Affinity purification IgG

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% Sodium azide,50% glycerol,pH7.3.

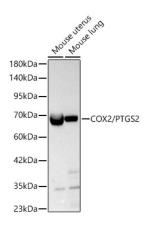

Western blot analysis of various lysates using COX2/PTGS2 Rabbit pAb (A1253) at 1:1000 dilution.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates / proteins: 25 µg per lane.


Blocking buffer: 3 % nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 90s.

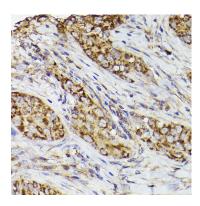

Immunohistochemistry analysis of paraffin-embedded Human breast cancer using COX2/PTGS2 Rabbit pAb (A1253) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

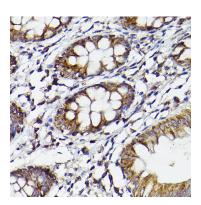

Immunohistochemistry analysis of paraffin-embedded Human colon using COX2/PTGS2 Rabbit pAb (A1253) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.