ABclonal www.abclonal.com

AKAP8 Rabbit pAb

Catalog No.: A17302

Basic Information

Observed MW

100kDa

Calculated MW

76kDa

Category

Polyclonal Antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human, Rat

Background

This gene encodes a member of the A-kinase anchor protein family. A-kinase anchor proteins are scaffold proteins that contain a binding domain for the RI/RII subunit of protein kinase A (PKA) and recruit PKA and other signaling molecules to specific subcellular locations. This gene encodes a nuclear A-kinase anchor protein that binds to the RII alpha subunit of PKA and may play a role in chromosome condensation during mitosis by targeting PKA and the condensin complex to chromatin. A pseudogene of this gene is located on the short arm of chromosome 9.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:100 - 1:200

Immunogen Information

Gene ID Swiss Prot 10270 O43823

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 592-692 of human AKAP8 (NP_005849.1)

Synonyms

AKAP-8; AKAP95; AKAP 95; AKAP-95; AKAP8

Contact

www.abclonal.com


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 50% glycerol, pH7.3.

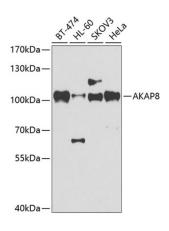

Western blot analysis of various lysates using AKAP8 Rabbit pAb (A17302) at 1:500 dilution.


Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 90s.

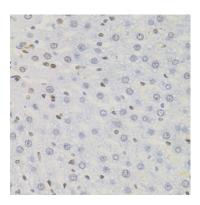

Immunohistochemistry analysis of paraffin-embedded rat liver using AKAP8 Rabbit pAb (A17302) (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

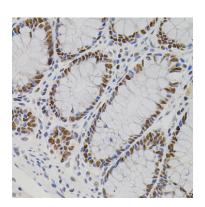

Immunohistochemistry analysis of paraffin-embedded human colon using AKAP8 Rabbit pAb (A17302) (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.