Leader in Biomolecular Solutions for Life Science

# Collagen X/COL10A1 Rabbit pAb

www.abclonal.com

ABclonal

Catalog No.: A18604 11 Publications

### **Basic Information**

### **Observed MW**

66kDa

### **Calculated MW**

66kDa

### Category

Polyclonal Antibody

### **Applications**

WB,IHC-P,ELISA

### **Cross-Reactivity**

Human, Rat

# Background

This gene encodes the alpha chain of type X collagen, a short chain collagen expressed by hypertrophic chondrocytes during endochondral ossification. Unlike type VIII collagen, the other short chain collagen, type X collagen is a homotrimer. Mutations in this gene are associated with Schmid type metaphyseal chondrodysplasia (SMCD) and Japanese type spondylometaphyseal dysplasia (SMD).

# **Recommended Dilutions**

WB 1:100 - 1:500

IHC-P 1:50 - 1:200

Recommended starting **ELISA** 

> concentration is 1 µg/mL. Please optimize the concentration based on your specific

assay requirements.

# **Immunogen Information**

**Gene ID Swiss Prot** 1300 Q03692

### **Immunogen**

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

# **Synonyms**

COL10A1; collagen alpha-1(X) chain; Collagen X/COL10A1

### Contact

www.abclonal.com €

### **Product Information**


Source **Isotype Purification** Rabbit Affinity purification IgG

#### Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% Sodium azide,50% glycerol,pH7.3.

# **Validation Data**

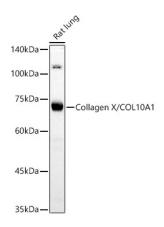

Western blot analysis of lysates from Rat lung, using Collagen X/COL10A1 Rabbit pAb (A18604) at 1:500 dilution.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 180s.

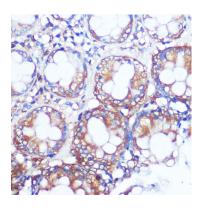

Immunohistochemistry analysis of paraffin-embedded Human colon using Collagen X/COL10A1 Rabbit pAb (A18604) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.