# SEC62 Rabbit mAb

Catalog No.: A24001 Recombinant 1 Publications



## **Basic Information**

#### **Observed MW**

50kDa

### **Calculated MW**

46kDa

### Category

SMab Recombinant Monoclonal Antibody

### **Applications**

WB,IHC-P,IF/ICC,ELISA

### **Cross-Reactivity**

Human, Mouse, Rat

### CloneNo number

ARC3193

# **Background**

The Sec61 complex is the central component of the protein translocation apparatus of the endoplasmic reticulum (ER) membrane. The protein encoded by this gene and SEC63 protein are found to be associated with ribosome-free SEC61 complex. It is speculated that Sec61-Sec62-Sec63 may perform post-translational protein translocation into the ER. The Sec61-Sec62-Sec63 complex might also perform the backward transport of ER proteins that are subject to the ubiquitin-proteasome-dependent degradation pathway. The encoded protein is an integral membrane protein located in the rough ER.

# **Recommended Dilutions**

**WB** 1:1000 - 1:2000

**IF/ICC** 1:200 - 1:400

IHC-P 1:200 - 1:800

**ELISA** Recommended starting

concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

### **Contact**

www.abclonal.com

## **Immunogen Information**

 Gene ID
 Swiss Prot

 7095
 Q99442

### **Immunogen**

Recombinant protein (or fragment). This information is considered to be commercially sensitive.

### **Synonyms**

HTP1; TP-1; Dtrp1; TLOC1; SEC62


### **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

#### Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

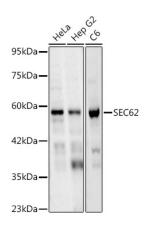

Western blot analysis of various lysates, using SEC62 Rabbit mAb (A24001) at 1:1000 dilution.

Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000


dilution

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 60s.

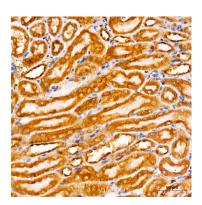


Immunohistochemistry analysis of paraffin-embedded Mouse kidney tissue using SEC62 Rabbit mAb (A24001) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

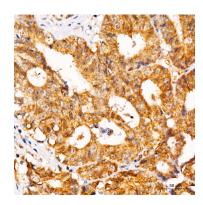


Immunohistochemistry analysis of paraffin-embedded Human colon carcinoma tissue using SEC62 Rabbit mAb (A24001) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

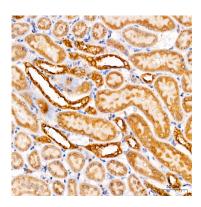


Immunohistochemistry analysis of paraffin-embedded Rat kidney tissue using SEC62 Rabbit mAb (A24001) at a dilution of 1:200 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IHC staining.

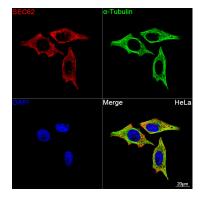

Confocal imaging of HeLa cells using SEC62 Rabbit mAb (A24001, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with  $\alpha$ -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG

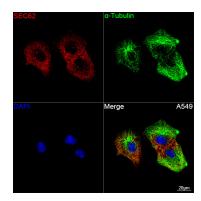

Confocal imaging of A549 cells using SEC62 Rabbit mAb (A24001, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with  $\alpha$ -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG

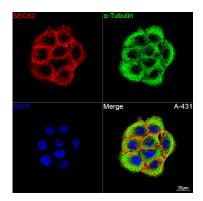

Confocal imaging of A-431 cells using SEC62 Rabbit mAb (A24001, dilution 1:200) followed by a further incubation with Cy3 Goat Anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). The cells were counterstained with  $\alpha$ -Tubulin Mouse mAb (AC012, dilution 1:400) followed by incubation with ABflo® 488-conjugated Goat Anti-Mouse IgG

## **Validation Data**

(H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x. (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x. (H+L) Ab (AS076, dilution 1:500) (Green). DAPI was used for nuclear staining (Blue). Objective: 100x.