# **VAMP8 Rabbit mAb**

Catalog No.: A4728 Recombinant



### **Basic Information**

#### **Observed MW**

15kDa

#### **Calculated MW**

11kDa

#### **Category**

SMab Recombinant Monoclonal Antibody

### **Applications**

WB,IHC-P,ELISA

#### **Cross-Reactivity**

Human, Mouse, Rat

#### CloneNo number

ARC1108

# **Background**

This gene encodes an integral membrane protein that belongs to the synaptobrevin/vesicle-associated membrane protein subfamily of soluble N-ethylmaleimide-sensitive factor attachment protein receptors (SNAREs). The encoded protein is involved in the fusion of synaptic vesicles with the presynaptic membrane.

# **Recommended Dilutions**

**WB** 1:500 - 1:1000

**IHC-P** 1:50 - 1:200

# **Immunogen Information**

**Gene ID**Swiss Prot

8673
Q9BV40

#### **Immunogen**

A synthetic peptide corresponding to a sequence within amino acids 1-100 of human VAMP8 (NP\_003752.2).

### **Synonyms**

EDB; VAMP-8; VAMP8

### **Contact**

www.abclonal.com

# **Product Information**


SourceIsotypePurificationRabbitIgGAffinity purification

# Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

# **Validation Data**

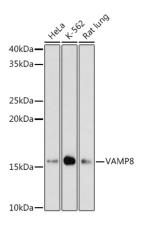

Western blot analysis of various lysates using VAMP8 Rabbit mAb (A4728) at 1:1000 dilution.


Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 180s.

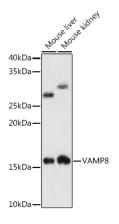


Western blot analysis of various lysates using VAMP8 Rabbit mAb (A4728) at 1:1000 dilution.

Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Enhanced Kit (RM00021).

Exposure time: 180s.

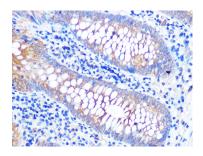

Immunohistochemistry analysis of paraffin-embedded Human appendix using VAMP8 Rabbit mAb (A4728) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.