

A5683

Leader in Biomolecular Solutions for Life Science



BNIP3 Rabbit pAb

Catalog No.: A5683

18 Publications

Basic Information

Observed MW

22-70kDa

Calculated MW

22kDa

Category

Polyclonal Antibody

Applications

WB,IHC-P,ELISA

Cross-Reactivity

Human,Mouse,Rat

Background

This gene encodes a mitochondrial protein that contains a BH3 domain and acts as a pro-apoptotic factor. The encoded protein interacts with anti-apoptotic proteins, including the E1B 19 kDa protein and Bcl2. This gene is silenced in tumors by DNA methylation.

Recommended Dilutions

WB 1:500 - 1:1000

IHC-P 1:50 - 1:200

ELISA Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID

664

Swiss Prot

Q12983

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 40-120 of human BNIP3 (NP_004043.3).

Synonyms

NIP3; HABON; BNIP3

Contact



www.abclonal.com

Product Information

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage


Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% Sodium azide, 50% glycerol, pH7.3.


Validation Data

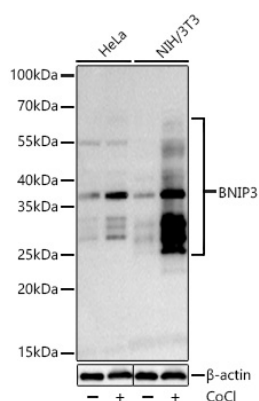

Western blot analysis of various lysates using BNIP3 Rabbit pAb (A5683) at 1:1000 dilution incubated overnight at 4°C. HeLa and NIH/3T3 cells were treated by Cobalt chloride (0.1 mM) at 37°C for 4 hours.
 Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.
 Lysates/proteins: 30 µg per lane.
 Blocking buffer: 3% nonfat dry milk in TBST.
 Detection: ECL Basic Kit (RM00020).
 Exposure time: 45s.

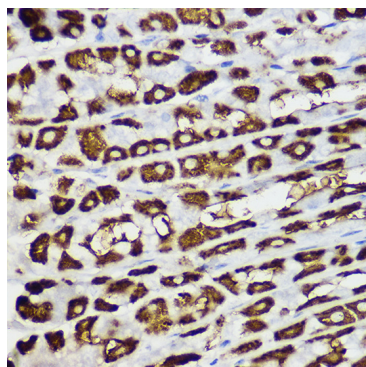

Immunohistochemistry analysis of paraffin-embedded Mouse stomach using BNIP3 Rabbit pAb (A5683) at dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

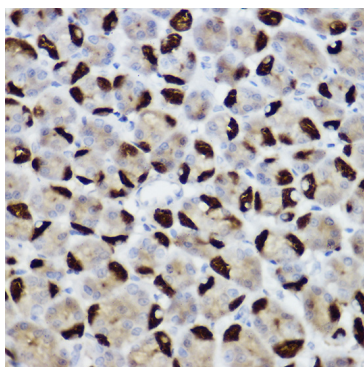

Immunohistochemistry analysis of paraffin-embedded Rat stomach using BNIP3 Rabbit pAb (A5683) at dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.