

RPT0002LQ

Leader in Biomolecular Solutions for Life Science



Recombinant *Schistosoma japonicum* GST-His Protein

Catalog No.: RPT0002LQ

Recombinant

Sequence Information

Species	Gene ID	Swiss Prot
<i>E. coli</i>		P08515

Tags

C-His

Synonyms

GST 26;Sj26 antigen;SjGST

Product Information

Source	Purification
<i>E. coli</i>	≥ 95 % as determined by SDS-PAGE.

Endotoxin

Please contact us for more information.

Formulation

Supplied as a 0.22 µm filtered solution in PBS, pH 7.4.

Reconstitution

Background

Genetic engineers have used glutathione S-transferase to create the GST gene fusion system. This system is used to purify and detect proteins of interest. In a GST gene fusion system, the GST sequence is incorporated into an expression vector alongside the gene sequence encoding the protein of interest. Induction of protein expression from the vector's promoter results in expression of a fusion protein: the protein of interest fused to the GST protein. This GST-fusion protein can then be purified from cells via its high affinity for glutathione. GST is commonly used to create fusion proteins. The tag has the size of 22 amino acids (roughly 26 kDa), which, compared to other tags like the Myc- or the FLAG-tag, is quite big. However, many commercially-available sources of GST-tagged plasmids include a thrombin domain for cleavage of the GST tag during protein purification.

Basic Information

Description

Recombinant *Schistosoma japonicum* GST-His Protein is produced by *E. coli* expression system. The target protein is expressed with sequence (Met1-Lys218) of *Schistosoma japonicum* GST-His fused with 6×His tag at the C-terminus.

Bio-Activity

Shipping

The product is shipped at ambient temperature. Upon receipt, store it immediately at the temperature recommended below.

Operational Notes

For your safety and health, please wear a lab coat and disposable gloves for handling.

Storage

Store at -70°C. This product is stable at ≤ -70°C for up to 1 year from the date of receipt. For optimal storage, aliquot into smaller quantities after centrifugation and store at recommended temperature. Avoid repeated freeze-thaw cycles. Avoid repeated freeze/thaw cycles.


Contact



www.abclonal.com

* For your safety and health, please wear a lab coat and disposable gloves when handling.

Validation Data

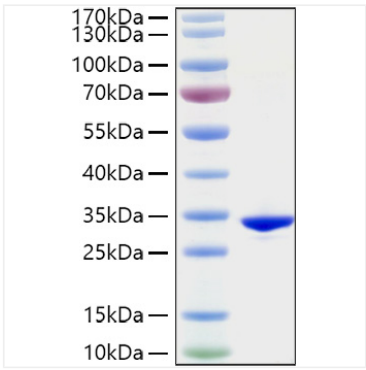

Recombinant *Schistosoma japonicum* GST-His Protein was determined by SDS-PAGE under reducing conditions with Coomassie Blue.